Ultherapy® Staggered Treatment Plan



Visible confidence is just a lift away

Start any way you want. Get a full Face, Neck & Décolletage Ultherapy® treatment in a single session, or break it up to suit your time and budget. Then, choose a dual-depth treatment or enhance with a third depth whenever you want.

1 FRESH EYES

TOTAL LINES = 100 +/- 20 LINES







• In patients with a thinner dermis, the frontalis bone may be more superficial than 4.5mm. In this case, you may forgo treatment with DS7-4.5 and add this to the number of lines delivered with the DS7-3.0

Standard Treatment 3.0 Depth 100 lines	
Enhance Your Treatment Add 1.5mm Depth 90 lines	
Total Price	



2 CONTOURED JAWLINE

TWO PHASES
IN TOTAL
430 lines

• Treating the submalar (lower cheek) and submentum has been shown to improve jawline definition and reduce submentum laxity¹

Standard Treatment 3.0 + 4.5mm Depth 430 lines	
Enhance Your Treatment Add 1.5mm Depth (enhance) 190 lines	
Total Price	

Submalar Treatment Total Lines: 430





Phase 2 (1 Week After)

Submentum Treatment Total Lines: 230






3

SMOOTH NECK LIFT

TOTAL LINES: 420

Standard Treatment 3.0 + 4.5mm Depth	
Enhance Your Treatment Add 1.5mm Depth	
Total Price	



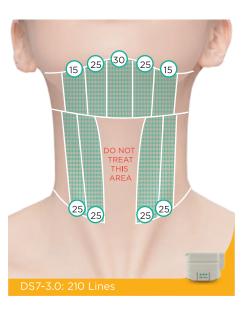

- Do not treat over the thyroid gland, thyroid cartilage or trachea
- Avoid treating directly into bone, as this may increase discomfort
- Avoid treating directly into major vessels

- Do not treat posterior to the preauricular line / tragus
- Pulls skin tight when sliding the transducer to prevent dragging the skin
- Refer to the IFU for additional treatment information

YOUTHFUL DÉCOLLETAGE

TOTAL LINES: 280

Standard Treatment 4.5 + 3.0 + 1.5mm Depth	
Total Price	





